

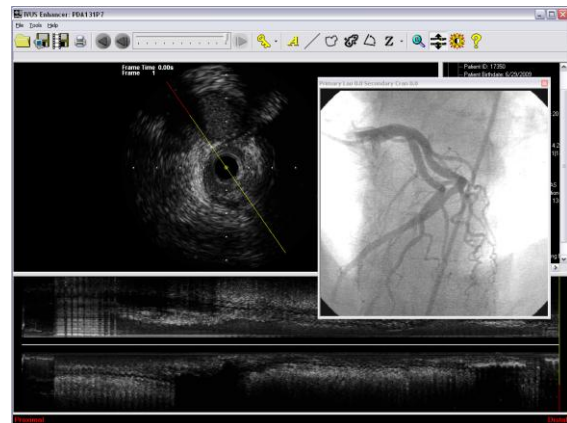
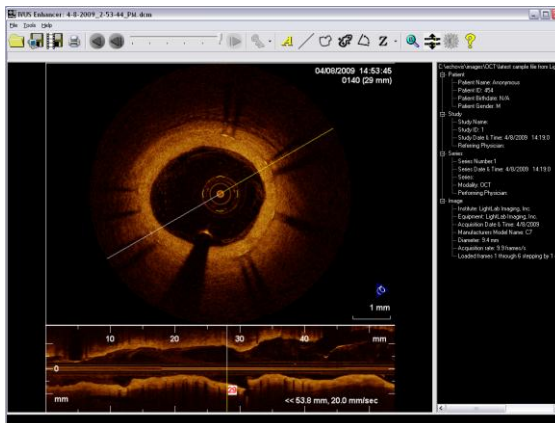


Enhancer: Multi-modality DICOM Review, Analysis, and Reporting

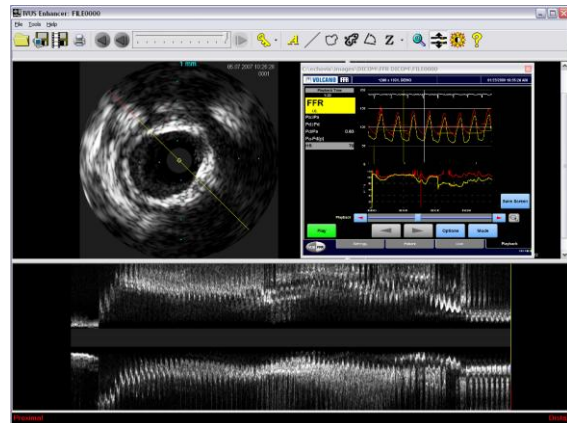

INDEC's **Enhancer** is a versatile, vendor-neutral software product for review and analysis of DICOM, echoPlaque, and Boston Scientific Galaxy™ files. **Enhancer** offers powerful capabilities including full length, real-time longitudinal display and rotation; image size, brightness, and contrast adjustment; two dimensional distance and area quantification; patient report generation; and the ability to export still images and video animations. Although primarily designed for use with IVUS images, **Enhancer** is compatible with most DICOM modalities including OCT (Optical Coherence Tomography), Ultrasound (US), and Angiography images. **Enhancer**'s Pic-in-Pic feature allows simultaneous IVUS and Angio viewing. Moreover, **Enhancer** is laptop compatible, so now image analysis can travel with you wherever you go.


Enhancer is compatible with most DICOM modalities including OCT (Optical Coherence Tomography), Ultrasound (US), and Angiography images. **Enhancer**'s Pic-in-Pic feature allows simultaneous IVUS and Angio viewing. Moreover, **Enhancer** is laptop compatible, so now image analysis can travel with you wherever you go.

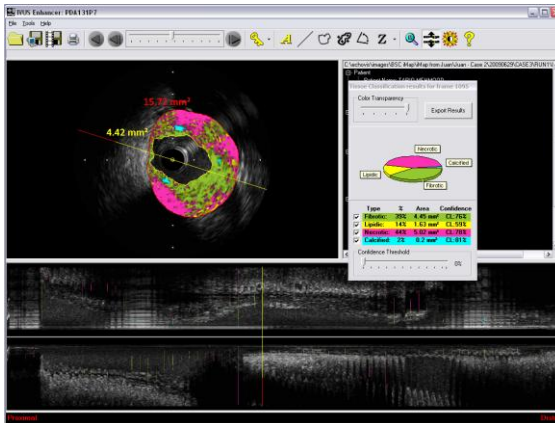


INDEC Systems, Inc. 4701 Patrick Henry Drive, Building 24, Suite 2401, Santa Clara, CA 95054

Tel: 1-408-986-1600 www.indecsystems.com

- Supports a variety of IVUS and OCT file formats including: DICOM, Boston Scientific Galaxy native files, and INDEC's proprietary echoPlaque files
- Supports wide range of DICOM modalities such as IVUS, US (Ultrasound), XA (Angiography), MR, and CT
- Simultaneous IVUS and Angio viewing with Pic-In-Pic
- Images automatically calibrated, patient demographic and study information automatically read
- Loads existing measurements, audio, ECG, and bookmarks
- Add new distance, area, and arc measurements or add annotations to the image
- Resize image or zoom in on specific regions
- Generate Patient Reports with a few simple clicks
- Export images and animations for presentations
- For use on desktop or laptop computers
- Fully-compatible with INDEC's Liberator database software - query and retrieve files from remote storage system

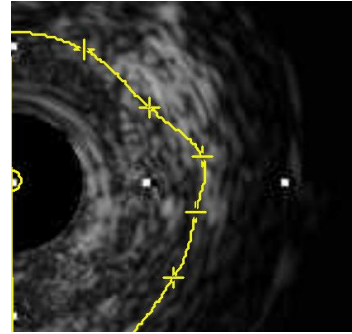

Versatile multi-modality image support

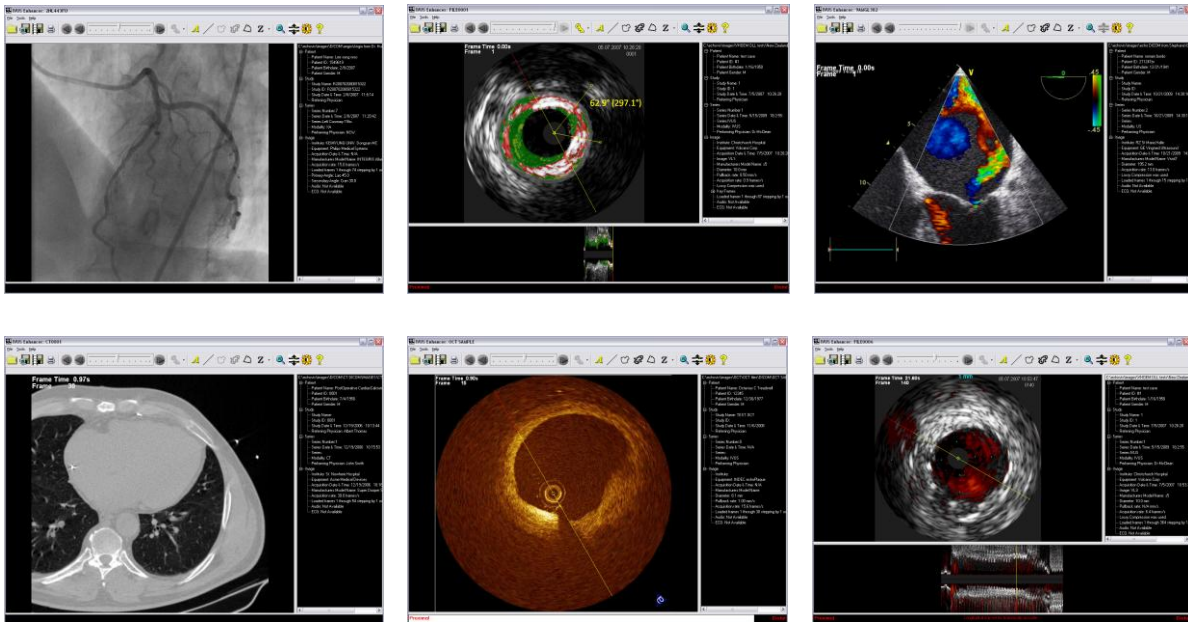

INDEC Systems, Inc. 4701 Patrick Henry Drive, Building 24, Suite 2401, Santa Clara, CA 95054

Tel: 1-408-986-1600 www.indecsystems.com